

# FLEX<sup>®</sup> Retractor

## Advanced Exposure for Advanced Procedures<sup>™</sup>



# FLEX<sup>®</sup> Retractor

The fast, easy way to confidently and consistently achieve advanced exposure

## KEY ADVANTAGES

- Uniquely designed for TORS, TLM and Flex<sup>®</sup> procedures
- Fully adjustable for optimal exposure in more patients
- Low profile, lightweight design easy to handle and position
- Rapidly swappable blades to avoid lengthy repositionings
- Comprehensive blade set for access to oropharynx, larynx and hypopharynx

Jaw spreader with paddle for secure tightening and quick release

Adjustable tongue blade angle enhances exposure and access options

Adjustable axial blade positioning to enhance tongue base exposure

Blades move independently from other retractor components and load and unload from top for easy and quick swaps and adjustments

Integrated smoke evacuation to efficiently maintain visualization


Adjustable cheek retraction for more flexible positioning

Adjustable tongue depressor insertion depth of 30mm increases versatility of blades and creates exceptionally low profile

Tooth brace options to accommodate varied dentition

Large working space opens up to 8cm and accommodates TORS, TLM and FLEX<sup>®</sup> procedures

Only the Flex<sup>®</sup> Retractor enables axial rotation of blade for unprecedented exposure of base of tongue

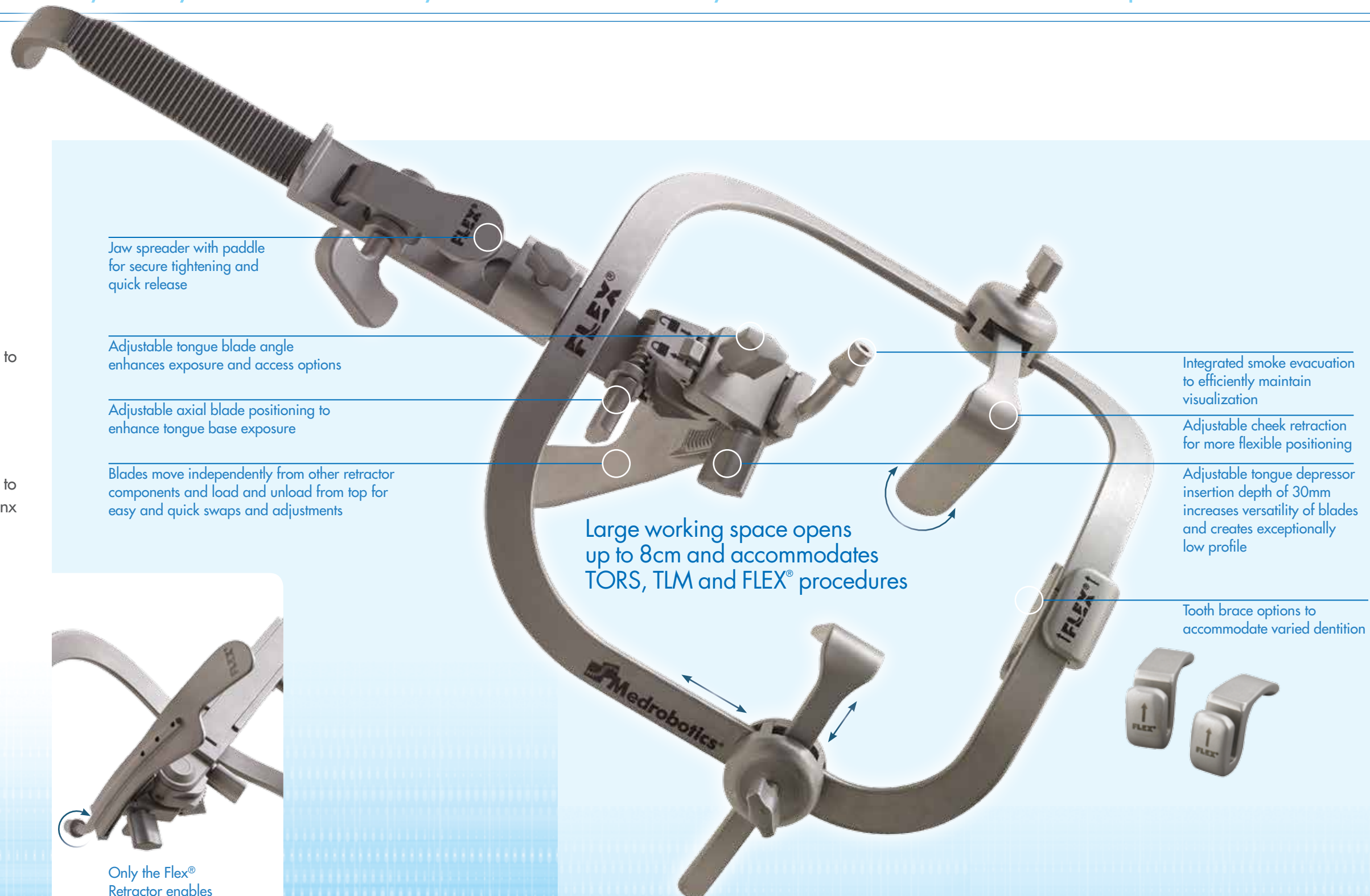

# More blades for more efficient procedures

Comprehensive array of easily swapped tongue blades facilitates the right exposure and procedural efficiency

## COMPLETE BLADE SET INCLUDED

### Tongue – Mandible Blades



Insertion Depth: 50-80mm  
Width: 32mm



Insertion Depth: 50-80mm  
Width: 38mm



Insertion Depth: 80-110mm  
Width: 32mm



Insertion Depth: 50-80mm  
Width: 40mm



Insertion Depth: 50-80mm  
Width: 40mm

### Laryngeal Blades



Insertion Depth: 90-120mm  
Width: 32mm



Insertion Depth: 90-120mm  
Width: 36mm